# PLANTA

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## AREA OF USE

---

With MF-20 massager you can enjoy massage at any time without leaving home. The device is used for massage of feet, ankles, calves, knees, thighs, hands. Massage with this massager has relaxing and stimulating effect, reduces tension and increases muscle tone, eliminates tension, helps to get rid of fatigue. The heating function stimulates blood circulation. MF-20 massager is ideal for people who stand on their feet for a long time at work and experience discomfort in this regard, as well as for those who lead a sedentary lifestyle.

This device is not a complete replacement for medical therapy or physiotherapy. If you are not sure that this massager is suitable for you, consult with your doctor.

The device is intended for personal use. It is not intended for commercial use of any kind, nor for use in medical institutions. To use all advantages of your massager carefully read the user manual and follow the rules of care and safety. This will extend the service life of your massager.

Save the user manual for future use or when transferring the device to another person. For more information visit our website: [www.planta.ru.com](http://www.planta.ru.com) .

**ATTENTION! THE MASSAGER IS NOT A MEDICAL PRODUCT!**

### Advantages of the MF-20 massager

- **Multifunctional** - MF-20 can be used for massage of different parts of body (feet, ankles, calves, thighs, knees, arms). Also this massager consists of two blocks - the upper wireless block and the bottom block. With the MF-20 you can massage several zones at the same time and also divide the blocks for two users.
- \* **3 Programs** - restorative, relaxing and intensive massage.
- \* **3 levels of Intensity** - you can adjust the strength of the Intensity yourself - low, medium or high - depending on your preferences.
- \* **Heat** - increase the effect of massage, improves blood flow and promotes muscle relaxation.
- \* **Massage rollers** – 3 different types of massage rollers effectively massage the entire foot from the fingertips to the heel affecting the acupuncture points.
- \* **Compression** - affects the massage area from all sides, creates 3D effect and increase the effect of massage rollers.

## SAFETY INSTRUCTIONS

---

- Before using of the massager read the user manual and safety instructions carefully.
- This device operates from power mains. That means that it is subject to the same safety rules as any other electrical appliance.
- Use the device only with the correct voltage in power mains.
- Ensure that all parts of the device are in good condition before use. If a defect or fault is revealed, turn off the device and immediately pull the plug out of the socket. Only specialists of authorized services can repair the device. Never open the device yourself. Separate components of this device cannot be disassembled, repaired or reinstalled.
- Cable must not touch sharp edges.
- Do not pull the cable or the device itself to disconnect or transport the device.
- Do not use the device more than 15minutes. If you do not use the device, disconnect it from the

power mains.

- This device is forbidden to be used by persons (especially children) with physical, sensory or mental disabilities, as well as persons who do not have sufficient experience or knowledge to handle the device, except when using the device under the control of the responsible person or after obtaining the necessary information.
- Ensure that children do not use the device as a toy.
- Never use the device in close proximity to heat sources.
- Use the device only in absolutely dry rooms, making sure that your body is also dry. It is forbidden to use the device in the bathroom, shower, above sink filled with water, in swimming pools, sauna and outdoors.
- Always unplug the device before cleaning. During cleaning make sure that liquids do not hit the switches and buttons. The device can only be used when it is completely dry.
- It is forbidden to disassemble the device. It has no parts that the user can repair.
- It is forbidden to remove the panels of the device, since the repair of internal components can only be carried out by qualified specialists.
- Do not insert needles or other sharp objects into the massager.
- When this device is disposed of make it unusable by cutting off the cable.
  
- Never use this device in conjunction with electro-medical equipment, electronic life support systems such as ventilators or portable medical devices such as an electrocardiograph. This device can cause disturbances in the operation of the above-mentioned medical equipment.
  
- Do not use the device if there are open wounds, burns, hematomas, tumors, etc. on your body.
- Do not use the device when under the influence of alcohol.
- Stop massage if you have unpleasant or painful feelings.
  
- If you experience limb pain during or after use of the device, stop using immediately and see your doctor.
- If the power cord is damaged, it must be replaced in the service center or by a qualified person to avoid danger.
- The device has heating surfaces. Persons who are sensitive to the effects of heat should take into account that the useful life of the device is 15 minutes and the device can only be used under supervision.

## Warning

Consult your doctor before using the device if :

- You are pregnant;
- You use a pacemaker and others;
- You suffer from one of the listed diseases: influenza, fever, varicose veins, thrombosis, hepatitis, diabetes, inflammation, ischias, loss of intervertebral discs, etc.;
- You undergo rehabilitation after surgery.

## Notes

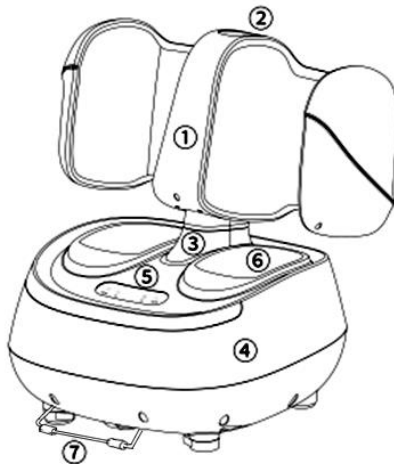
If you use the device just before going to sleep, remember that massage can have a stimulating effect and cause insomnia.

## DESCRIPTION OF THE DEVICE

---


### Box contents:

1. Massager which consists of 2 blocks – main foot part and removable wireless block
2. User Manual with warranty card

1. Upper removable wireless block
2. Control panel of Upper removable wireless block
3. Shaft for connecting of Upper block to Bottom block
4. Bottom foot block
5. Control panel of Bottom foot block
6. Place for foot placement inside the massager
7. Base column support

## CONTROL PANEL OF UPPER WIRELESS BLOCK

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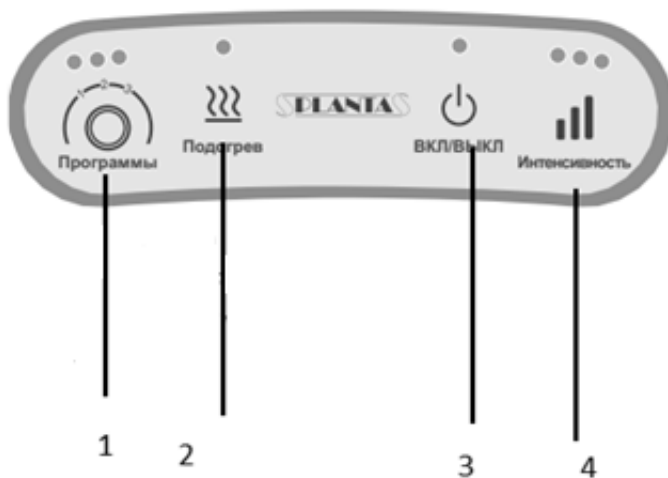Using this button you can select one of 3 preset programs. Program 1 starts automatically. Pressing the program button switches the work of the massager between 3 programs - which program is selected is shown by the light indication above the button.

## Description of programs for the bottom foot block

### Program 1

Restorative massage

Multiple massage of feet by rollers. During this program air compression in upper and bottom foot blocks helps to relax the muscles and increase blood flow.

### Program 2

Relaxing massage

Massage by air compression

### Program 3

Intensive massage

Compression airbags fix the body area, this compression with combination of the rollers has a profound effect on the active points of the feet. During this program special attention is paid to the tendons.

2. Heat on and off button (Подогрев) 

Heat is automatically turned on when the massager is turned on. Pressing the button turn off the heating function. The indicator light above the button lights up when the heating function is activated.

3. On and off button (ВКЛ/ВЫКЛ) 

Use on and off button to turn on or off the massager. Pressing this button during operation stops the program cycle and turn off massager.

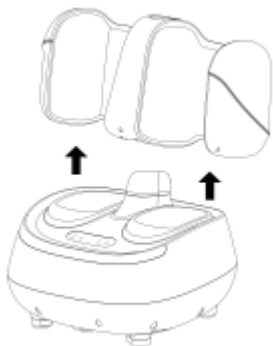
4. Button for choosing of Intensity level (Интенсивность) 

You can choose one of 3 intensity levels. When turned on the massager starts working at medium 2nd intensity level. Pressing this button once increases the intensity to the 3rd level. Further pressing will return the intensity of the massage to the 1st level. The selected intensity level is displayed by a light indicator above the button.

## Upper removable wireless block

### Separation of the upper wireless block

Hold the upper block with 2 hands and lift it up to separate the upper block from the bottom block

### Changing position of the upper wireless block

Rotate the upper block clockwise or counterclockwise to easily adjust the position of upper block

**Attention! During separation, connection or changing the position of the upper wireless block the massager must be turned off.**

## OPERATING THE APPLIANCE

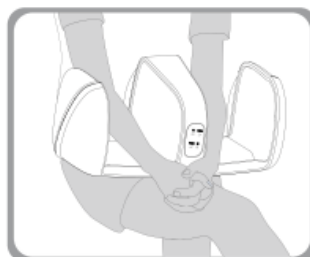
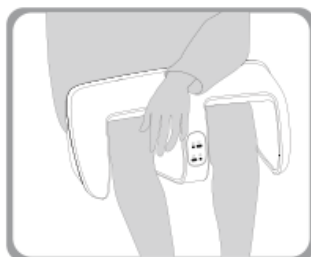
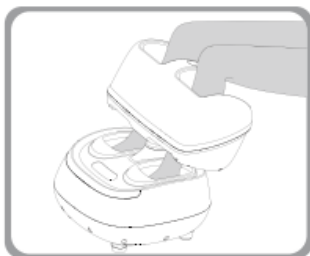
---

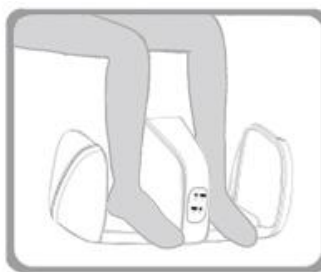
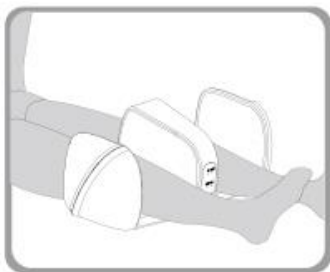
**Attention!** With this device you can massage several zones at the same time and also you can divide the blocks for two users. When using the upper and bottom blocks together you can control the massager from any control panel, when using the blocks separately control from each panel separately.

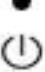
- Inspect the massager. Make sure that there is no damage of the massager and the power cord.
- Connect the device to a 220 V power mains. When you use only the wireless removable block separately, it is possible not to connect the bottom part to the power mains.
- Take a comfortable position. Place your feet inside the bottom block of the massager or place the removable block on the selected part of the body for massage. If necessary adjust the angle of inclination of the bottom foot block using the footrest.

- Choose one of the following options for using the massager.

• **ON/OFF.** Press button ВКЛ/Выкл . The indicator  will light up on the control panel, the massager will start working at medium intensity mode.

To turn off the massager, press button ВКЛ/Выкл , the operation indicator  will go out.

#### • CHOOSING OF INTENSITY LEVEL.

When turned on the massager starts working at medium intensity level.

Press button Интенсивность  to control the intensity of the air pressure and the speed of rotation of the rollers. You can choose your preferred level of intensity (it can be light foot treatment, medium massage or strong when the muscles are very tired). It is recommended to try the low level first and gradually increase the power.

When you press the button Интенсивность , the medium intensity switches to high intensity, the next press sets the low intensity, then the medium again and so on. The corresponding indicator on the control panel lights up.

#### • PROGRAMS.

When turned on, the massager will start working according to the Program 1.

Press button Программы  to control the programs. Every time you press the Program Button, there is a switch between Program 1 -> Program 2-> Program 3 -> Program 1 and so on, and the corresponding indicator lights up.

#### NOTES:

•The device switches off automatically after 15 minutes for safety reasons. This function should not be

considered as a replacement for the turn off function. Do not forget to disconnect the appliance from the mains if you do not use it.

- Do not try to stand on this massager. Use the bottom block as well as the removable wireless unit only when you are sitting. The removable wireless block can be used when you are lying down.
- Do not use the massage or heat functions continuously for more than 15 minutes. The beneficial effect of the massage is 15 minutes - this corresponds to the optimal duration of the massage which has a beneficial effect on the body.

### **Charging of Upper removable wireless block**

1. Connect the upper and bottom massager block, plug the massager into the power mains;
2. When charging the massager the power indicator flashes only on the upper block. When the upper block is fully charged, the power indicator turns off.
3. It takes about 4 hours to fully charge the device.
4. Do not use the massager while charging the upper wireless block.

### **Attention!**

- If defects are revealed in the operation of the device, immediately stop using it and contact the seller
- For safety reasons, all repair and maintenance work on the device must be carried out only by authorized representatives of the manufacturer. Do not try to fix the problem yourself!

## **CARE AND STORAGE**

---

- Turn off the device and let it cool down before cleaning.
- Use cloth wetted with soft detergent for cleaning.

**Attention! Never use abrasive or aggressive cleaning agents and do not put the device into the water.**

- The device can be used again after it is completely dry.

**Attention! It is forbidden to remove the panels of the device or open the fasteners because the repair of internal components can only be carried out by qualified specialists.**

- Removable covers on the upper and bottom block can be unfastened with zippers and then washed by hands .Let the covers completely dry after washing before installing them back on the appliance.
- If you do not use the massager, store it in a dry place away from direct sunlight.
- Let device to cool completely before folding it.
- Do not place heavy objects on the device during storage otherwise parts of the device may break.
- Do not allow the device to fall. This may cause damage.

## POSSIBLE PROBLEMS AND HOW TO FIX THEM

Problem	Possible cause	How to fix
Device does not switch on.	The device is not connected to a power source.	Make sure that the device is connected to power.
	The device was switched on not correctly.	Press on and off button ВКЛ/ВЫКЛ
The device suddenly stopped working.	You use the massager more than 15 minutes.	It is normal. The device works properly.
Upper wireless block is stopped working/does not turn on	Upper wireless block is discharged	Charge upper wireless block
The device works not correctly.	The device was kept at low temperature for some time.	Leave the massager for 2 hours at room temperature, then use it.
The sound of escaping air is heard.	There is correction of the pressure in the air chambers.	It is normal situation for working of the massager.

### Disposal of electrical equipment

In the interests of the environment, after the end of the service life, the device cannot be disposed of together with household waste. Disposal shall be carried out through appropriate collection and processing points for electrical and electronic equipment. If you have questions, contact the local municipal service responsible for waste disposal.

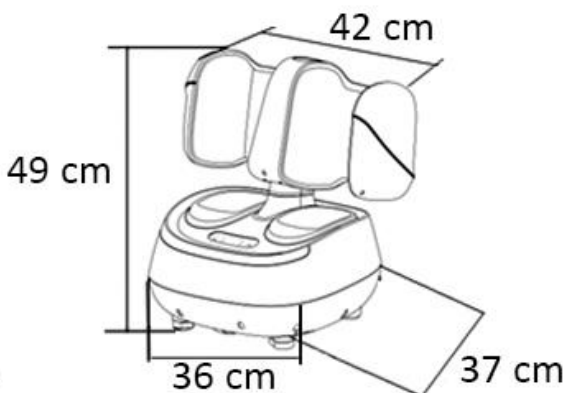
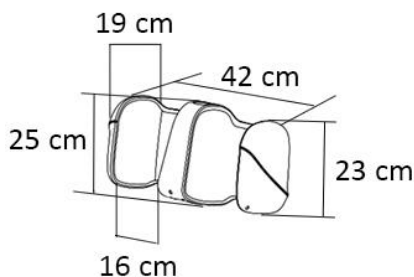
## TECHNICAL SPECIFICATION

Model: MF-20

Rated voltage: 220V / 50Hz

Power: 65W


Dimensions:

Weight: 8,9 kg (3,4 kg upper wireless block and 5,5 kg bottom block)

Power cord size: 1.83 meters

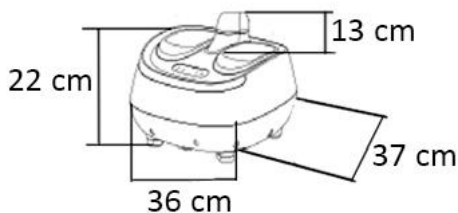Device material: ABS plastic, polyester

	High intensity	Medium intensity	Law intensity
Bottom block	76 kPa (0,57 mm Hg)	74 kPa (0,56 mm Hg)	71 kPa (0,53 mm Hg)
Upper block	73 kPa (0,55 mm Hg)	71 kPa (0,53 mm Hg)	69 kPa (0,52 mm Hg)

The manufacturer reserves the right to change parts of the product without prior notice for different schedule lines without affecting the basic technical parameters of the product. This can entail a change in the weight and dimensions of the product, but not more than +/- 5-10%.

# WARRANTY CARD

<p><b>Warranty Terms</b></p> <p>1.If a fault is revealed during the Warranty Period, the Product shall be exchanged with the retailer that sold the Product.</p> <p>2. Warranty replacement shall be carried out if the warranty card is filled in and subject to the operating rules described in the instruction manual.</p> <p>3.Product warranty does not apply in the following cases:</p> <ul style="list-style-type: none"><li>- mechanical damages;</li><li>- failure of the article due to ingress of foreign objects and liquids, insects, etc.;</li><li>- use of the article in conditions and modes other than domestic.</li></ul> <p>4.The warranty does not apply to accessories and components.</p> <p>5. The warranty also loses force if during the warranty period the repair of the faulty product was carried out by non-authorized persons.</p> <p>6.The warranty period is 1 year.</p> <p>7.The service life is 3 years.</p> <p>The date of production is indicated on the package.</p> <p>Manufacturer: XIAMEN COMFORT SCIENCE &amp; TECHNOLOGY GROUP CO., LTD</p> <p>Address: No. 168 QianPu Road,Siming Zone, Xiamen 361009, P.R.China</p> <p>Importer/Authorized Organization: LLC "PALLADA" 111675,Russia, Moscow, Rudnevka Street, House 4, P IV K 1.2 OF 1 RM 2/5, Russia</p> <p>Customer service support: <a href="https://www.planta.ru.com">https://www.planta.ru.com</a></p> <p>+7 (495) 137-80-32</p> <p>Product is certified</p> 	<p>WARRANTY CARD</p> <p>Product: _____</p> <p>Model: _____</p> <p>Date of sale: _____</p> <p>Selling company information: _____</p> <p>Product received, claims against configuration and appearance I do not mean, I agree with the terms of the guarantee. F.I.O. and buyer's signature: _____</p> <p>Opened the package, checked it and sold it. F.I.O. of the seller: _____</p> <p>Stamp</p>
---	--